

**Joint ICTP-IAEA Workshop on Monte Carlo Radiation Transport  
and Associated Data Needs for Medical Applications**

28 October – 8 November 2024

ICTP, Trieste, Italy

## Lecture 23

---

# egs++ particle sources

**Reid Townson**

Metrology Research Centre  
National Research Council Canada

---



Government  
of Canada

Gouvernement  
du Canada



# Particle sources

Particle sources within a Monte Carlo code provide the following initial properties for particles:

- position  $\vec{x}$
- direction  $\vec{u}$
- energy  $E$
- charge  $q$
- statistical weight  $w$
- latch variable
- time index

Some of these variables must be sampled according to a probability distribution, e.g., to reproduce the energy distribution of incident particles or their position in space for an extended source.

# The egs++ particle sources

- EGSnrc C++ source classes deliver the same functionality as all other EGSnrc source functions (e.g., sources for RZ applications), but offer much more flexibility.
- Any egs++ source class provides a central function `getNextParticle()`
  - samples particle position, direction,  $E$ ,  $q$  etc.
  - returns the number of independent particles generated so far
- The `getNextParticle()` function is called within every egs++ application

```
int EGS_Application::simulateSingleShower() {  
    ...  
    current_case =  
        source->getNextParticle(rndm,p.q,p.latch,p.E,p.wt,p.x,p.u);  
    ... }  

```

# Common source input syntax

- The source definition is placed in the `.egsinp` file, between delimiters:

```
:start source definition:
```

```
...
```

```
:stop source definition:
```

- A source object is defined using

```
:start source:
```

```
    name      = name_of_source
```


```
    library = source_library_name
```

```
    ( other source-specific input ...)
```

```
:stop source:
```

- The final simulation source can be made up of several sources
- The source definition section must contain the following key to select the simulation source: `simulation source = name_of_source`

# The EGS\_BaseSimpleSource class

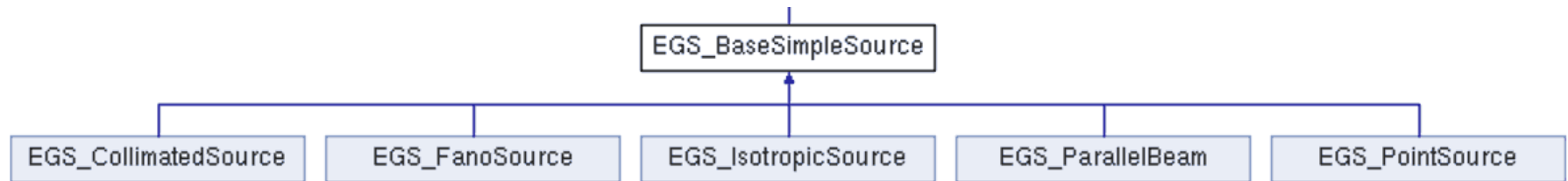

- Simple sources share the following properties:
  - deliver a fixed charge,
  - sample energy from a distribution using a spectrum object
  - use method `getPositionDirection()`: delivering  $\vec{x}, \vec{u}, w$
- Different types (parallel beam, collimated source, etc.) are derived from `EGS_BaseSimpleSource` and implement own `getPositionDirection()` functions  $\Rightarrow$  definition in input files is source-specific

**Source energy spectrum**

Source shape

Source type

# Source energy definition: Basics

Definition of the energy distribution is common to all and must be contained within

```
:start spectrum:
```

```
.
```

```
.
```

```
:stop spectrum:
```

The particle's energy can be either monoenergetic or sampled from an energy distribution (spectrum).

- Monoenergetic spectrum with a single discrete energy

```
:start spectrum:
```

```
  type = monoenergetic
```

```
  energy = the kinetic energy in MeV
```

```
:stop spectrum:
```

# Source energy definition: polyenergetic sources

- Uniform distribution within defined energy boundaries

```
type = uniform  
range = minimum and maximum energy  
or  
minimum energy = Emin  
maximum energy = Emax
```

- Tabulated spectrum: energies and probabilities in a file (e.g. \$HEN HOUSE/spectra)

```
type = tabulated spectrum  
spectrum file = absolute path to a spectrum file
```

or inline with

```
type = tabulated spectrum  
energies = list of discrete energies  
           or bin edges  
probabilities = list of probabilities  
spectrum type = 0 or 1 or 2 or 3
```


Note: probabilities do not need to be normalized!



# Tabulated spectrum types

`spectrum type` (or `MODE` within the file) defines one of four possible interpretations:

- Histogram spectrum (type = 0,1): series of energy bins with different uniform probabilities, considered as cts per bin (type = 0) or cts per MeV (type = 1)
- Line spectrum (type = 2): series of discrete energies with corresponding probabilities
- Interpolated spectrum (type = 3): probabilities considered to be at bin edges and linear interpolation between bin edges

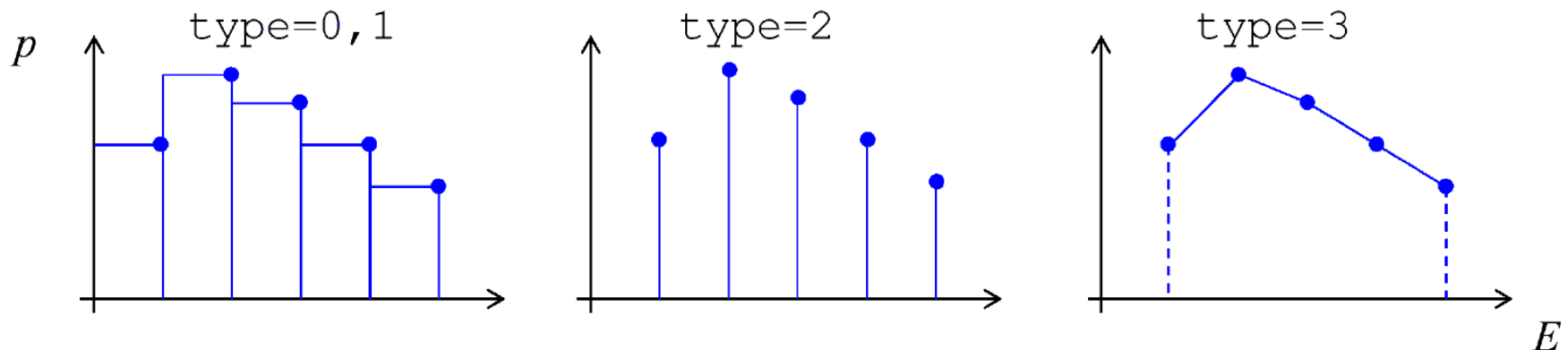

# Source energy definition: Gaussian sources

- Gaussian spectrum with energies sampled from a normal distribution

:start spectrum:

type = Gaussian

energy = the mean kinetic energy in MeV


sigma = the sigma of the spectrum

or

fwhm = the full-width-at-half-maximum of the spectrum

:stop spectrum:

- Double Gaussian spectrum: with two values for  $\sigma$ /fwhm, defining shape for energies less/larger than mean, joint smoothly at the mean energy.

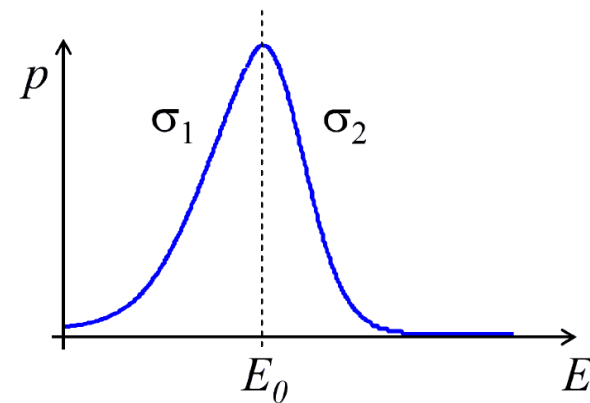


Source energy spectrum

**Source shape**

Source type

# EGS\_Shapes

- Sources derived from `EGS_BaseSimpleSource` require some method for the determination of direction and position of a particle:
  - position  $\vec{x}_0$  is fixed or picked randomly within a certain region in space
  - direction  $\vec{u}$  can be isotropic around  $\vec{x}_0$ , have fixed direction or aim to another point  $\vec{x}$ , thus  $\vec{u} = (\vec{x} - \vec{x}_0) / |\vec{x} - \vec{x}_0|$
- `EGS_shapes`
  - deliver random positions in space
  - internal function `getPoint()` samples and returns  $\vec{x}$
  - surface shapes deliver  $\vec{u}$  with the `getPointSourceDirection()` function, called with  $\vec{x}_0$

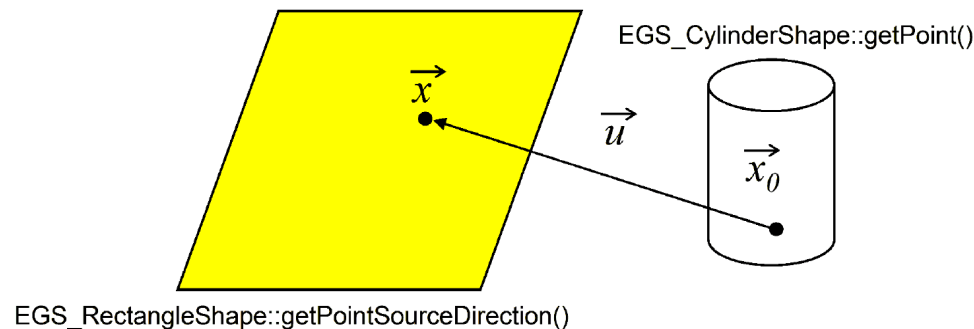

# EGS\_Shapes types

class	<b>EGS_BaseShape</b>	Base shape class. All shapes in the EGSnrc C++ class library are derived from <b>EGS_BaseShape</b> . <a href="#">More...</a>
class	<b>EGS_SurfaceShape</b>	A surface shape. <a href="#">More...</a>
class	<b>EGS_PointShape</b>	A point shape. This is the simplest shape possible: it simply always returns the same point. <a href="#">More...</a>
class	<b>EGS_BoxShape</b>	A box shape. <a href="#">More...</a>
class	<b>EGS_SphereShape</b>	A sphere shape. <a href="#">More...</a>
class	<b>EGS_CylinderShape</b>	A cylinder shape. <a href="#">More...</a>
class	<b>EGS_CircleShape</b>	A circle shape. <a href="#">More...</a>
class	<b>EGS_EllipseShape</b>	An elliptical shape. <a href="#">More...</a>
class	<b>EGS_ExtendedShape</b>	An extended shape. <a href="#">More...</a>
class	<b>EGS_GaussianShape</b>	A Gaussian shape. <a href="#">More...</a>
class	<b>EGS_LineShape</b>	A line shape. <a href="#">More...</a>
class	<b>EGS_TriangleShape</b>	A triangular shape. <a href="#">More...</a>
class	<b>EGS_PolygonShape</b>	A polygon shape. <a href="#">More...</a>
class	<b>EGS_RectangleShape</b>	A rectangular shape. <a href="#">More...</a>
class	<b>EGS_RectangularRing</b>	A "rectangular ring". <a href="#">More...</a>
class	<b>EGS_ShapeCollection</b>	A shape collection. <a href="#">More...</a>
class	<b>EGS_VoxelizedShape</b>	A "voxelized shape". <a href="#">More...</a>

# EGS\_Shapes specifics

- All have their own specific key values (see PIRS898).
- Some shapes are available as dso/dll and require the library key, others are directly compiled into egsp and defined with the type key.

- examples:

```
:start shape:  
    type = point  
    position = px, py, pz  
:stop shape:
```

or

```
:start shape:  
    library = egs_rectangle  
    rectangle = x1 y1 x2 y2 # left-upper, right-lower corners  
:stop shape:
```

- All shapes are placed at the origin by default and might be translated/rotated with the proper transformation.

## **EGS\_Shapes are not physical materials**

- The shapes just define the position sampling space for sources and do not have any physical manifestation in the geometry
- Notice: shapes do not have a material assigned!
- Source shapes can overlap geometries and are completely independent
- They are just a surface or volume within which source positions can be easily sampled

Source energy spectrum

Source shape

**Source type**



# EGS\_PointSource

The simplest type of source ...

```
:start source:
  library = egs_point_source
  name = some_name
  position = Px, Py, Pz
  :start spectrum:
    definition of the spectrum
  :stop spectrum:
  charge = -1 or 0 or 1 for electrons or photons or positrons
:stop source:
```

# EGS\_PointSource: Example

```
:start source definition:
  :start source:
    library      = egs_point_source
    name         = my_source
    position     = 0 0 0
    charge       = 0

    :start spectrum:
      type       = monoenergetic
      energy     = 1
    :stop spectrum:
  :stop source:

  simulation source = my_source

:stop source definition:
```

# EGS\_ParallelBeam

- `EGS_ParallelBeam` delivers particles, all having same direction with random positions within *any* shape

- Usually planar shapes are used, but 3D-shapes are possible.

`:start source:`

`library = egs_parallel_beam`

`name = some_name`

`:start shape:`

`definition of the shape`

`:stop shape:`

`:start spectrum:`


`definition of the spectrum`

`:stop spectrum:`

`direction = Ux Uy Uz`

`charge = -1 or 0 or 1 for electrons or photons or positrons`


`:stop source:`

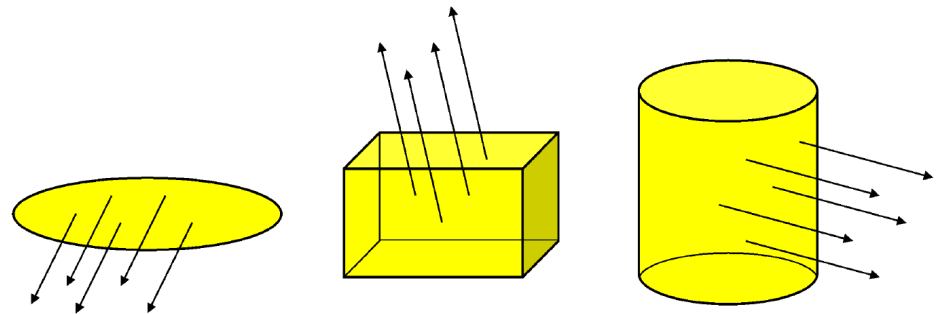

- Functionality of SOURCE NUMBER=0,2,10,13 from the RZ applications and isource=0,1 in `DOSXYZnrc` can be reproduced.

# EGS\_IsotropicSource

- [EGS\\_IsotropicSource](#) delivers particles with directions uniformly distributed in  $4\pi$  emitted from any shape.
- In case of 3D-shape, particles are emitted from any point within shape.
- Further, you can use a previously defined egsp geometry (or some selected regions of it) for sampling particle positions
  - IncludeAll, ExcludeAll, IncludeSelected or ExcludeSelected  
⇒ very complex isotropic sources shapes can be realized!
- Remember, the source shapes do not exist in the geometry and have no material.
- Source can reproduce the functionality of the SOURCE NUMBER=3 of the RZ-codes and isource=6 of DOSXYZnrc.

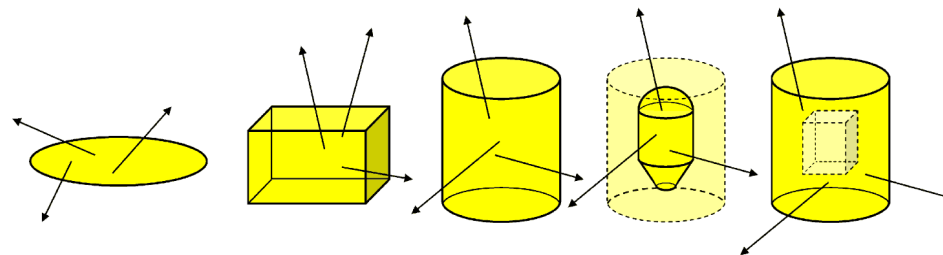


## EGS\_IsotropicSource: input definition

```
:start source:
  library = egs_isotropic_source
  name = some_name
  :start shape:
    definition of the shape
  :stop shape:
  :start spectrum:
    definition of the spectrum
  :stop spectrum:
  geometry = some_geometry
  region selection = IncludeAll or ExcludeAll or
                    IncludeSelected or ExcludeSelected
  #selected regions = selected regions to exclude or include
  charge = -1 or 0 or 1 #for electrons or photons or positrons
:stop source:
```

# EGS\_CollimatedSource

- `EGS_CollimatedSource` is an isotropic source, collimated to irradiate only a certain area or solid angle.

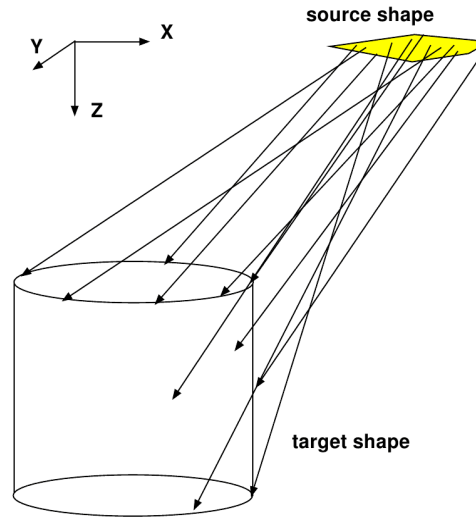


- Internally, `EGS_CollimatedSource` calls `getPoint()` function of a “source shape” and `getPointSourceDirection()` of a “target shape”  $\Rightarrow$  only shapes supporting this function can be used as targets.

## EGS\_CollimatedSource: example

```
:start source:
  name          = the_zero_collimated_source
  library        = egs_collimated_source
  distance       = 100
  charge         = 0
:start source shape:
  type          = point
  position       = 0, 0, -100
:stop source shape:
:start target shape:
  library        = egs_rectangle
  rectangle      = -15 -15 15 15
  # (in z=0 plane; use an affine transformation to change)
:stop target shape:
:start spectrum:
  type          = monoenergetic
  energy         = 0.06
:stop spectrum:
:stop source:
```

# EGS\_BeamSource

- `EGS_BeamSource` uses a `BEAMnrc` application, pre-compiled into shared library (dso/dll).
- Particles are extracted and used as source particles upon crossing a defined scoring-plane within `BEAMnrc`.
- When starting the simulation a complete `BEAMnrc` simulation is initialized.
- Internally, a container provides single particles of a primary history. When the container is empty, a new shower within `BEAMnrc` is called, filling the container.
- `getNextParticle()` function of `EGS_BeamSource` takes single particles from container.

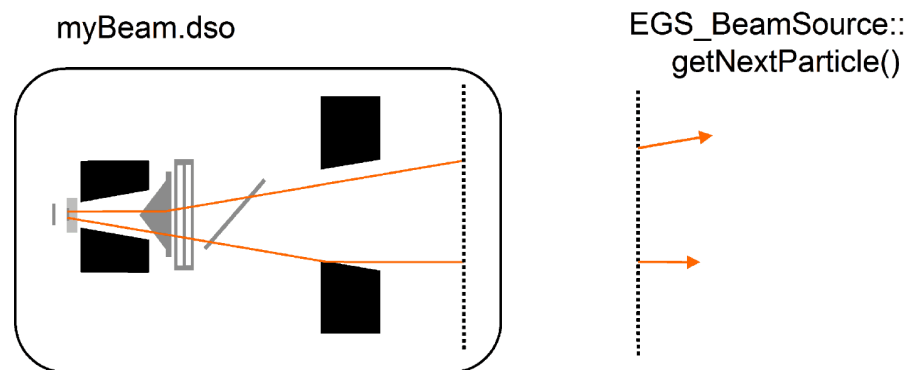




# EGS\_BeamSource: definition input

```
:start source:  
  library = egs_beam_source  
  beam code = BEAMnrc application name  
  pegs file = PEGS file name for BEAMnrc simulation  
  input file = BEAMnrc input file name  
  cutout = x1 x2 y1 y2 (optional)  
  particle type = 'all' or 'electrons' or 'photons'  
  or 'positrons' or 'charged' (optional)  
  weight window = wmin wmax (optional) # use to reject "fat" particles  
:stop source:
```

**BEWARE:** When restarting calculations using a BEAMnrc shared library source, make sure that the [restart](#) calculation option is defined in both the source and the application input files.

# EGS\_PhaseSpaceSource

A phase-space file source reads and delivers particles from a BEAMnrc phase-space file. Because the phase-space file only contains the x- and y- positions, the z-position is set to 0.

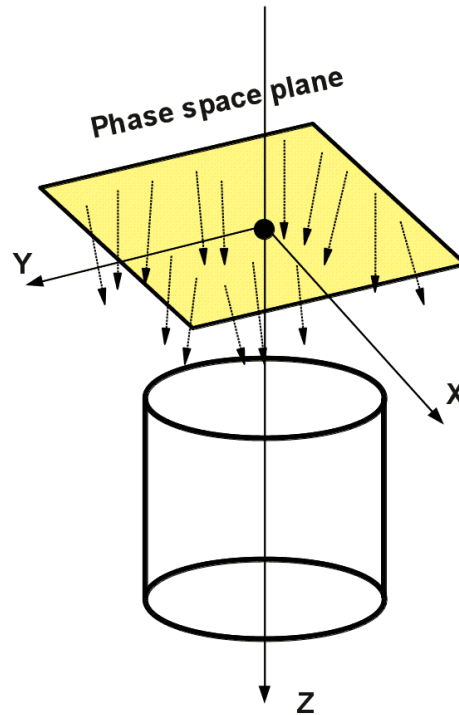

# EGS\_PhaseSpaceSource

A phase-space file source is defined as follows:

```
:start source:  
  library = egs_phsp_source  
  name = some_name  
  phase space file = phsp file name  
  particle type = photons # or electrons, positrons, all, or charged  
  cutout = x1 y1 x2 y2 (optional)  
  reuse photons = number of times to RECYCLE a photon  
  reuse electrons = number of times to RECYCLE a charged particle  
:stop source:
```

# EGS\_TransformedSource

A transformed source takes a particle from any other source and applies an **affine transformation** to the position and the rotation of to the direction.

A transformed source is defined as follows:


```
:start source:
  library = egs_transformed source
  name = some_name
  source name = the name of a previously defined source
  :start transformation:
    translation = tx, ty, tz
    rotation = 2, 3 or 9 floating point numbers
    or
    rotation vector = 3 floating point numbers
  :stop transformation:
:stop source:
```

## EGS\_PhaseSpaceSource redux

- Note that a phase-space source can be used as the source in a transformed source permitting in this way arbitrary transformations to be applied to the particle positions and directions.
- It is also worth noting that, together with a transformation, the phase-space source can reproduce the functionality of any phase-space file based source in the RZ series of applications and in DOSXYZnrc.


# EGS\_RadionuclideSource

- A new source of particles, following the decay scheme of a radionuclide.
- `EGS_RadionuclideSource` extends `EGS_BaseSource`, and allows for the correlated emissions of photons, electrons and positrons.

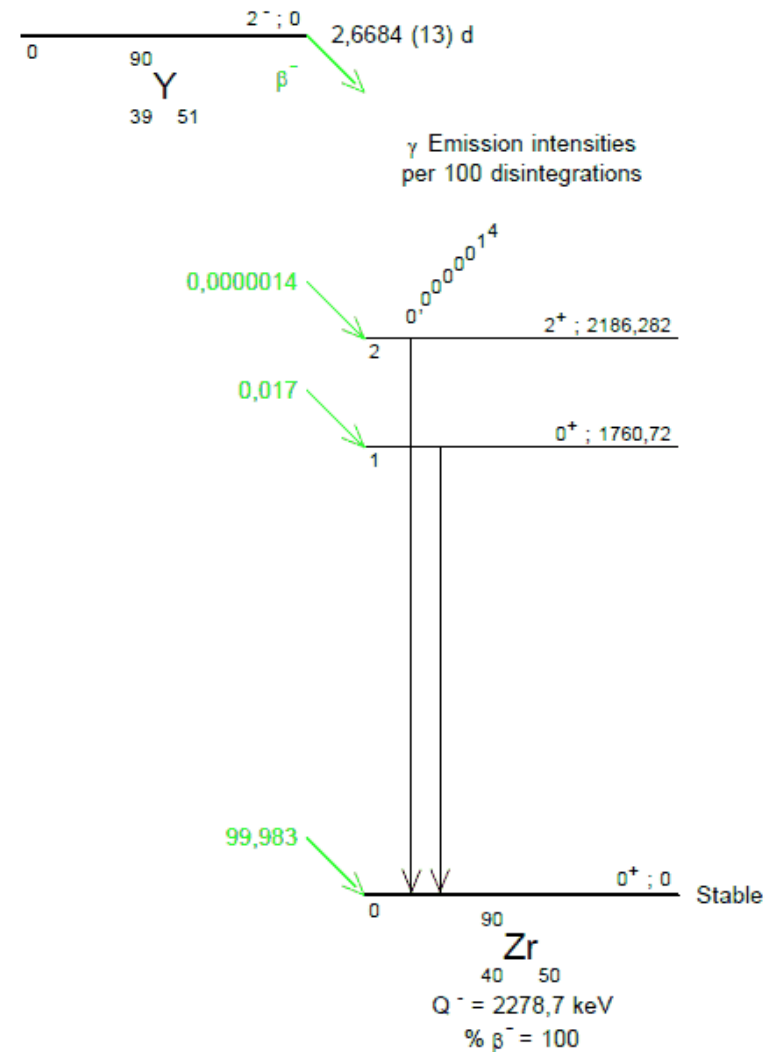

Nuclear data from LNHB: [http://www.nucleide.org/DDEP\\_WG/DDEPdata.htm](http://www.nucleide.org/DDEP_WG/DDEPdata.htm)

# EGS\_RadionuclideSource (cont.)

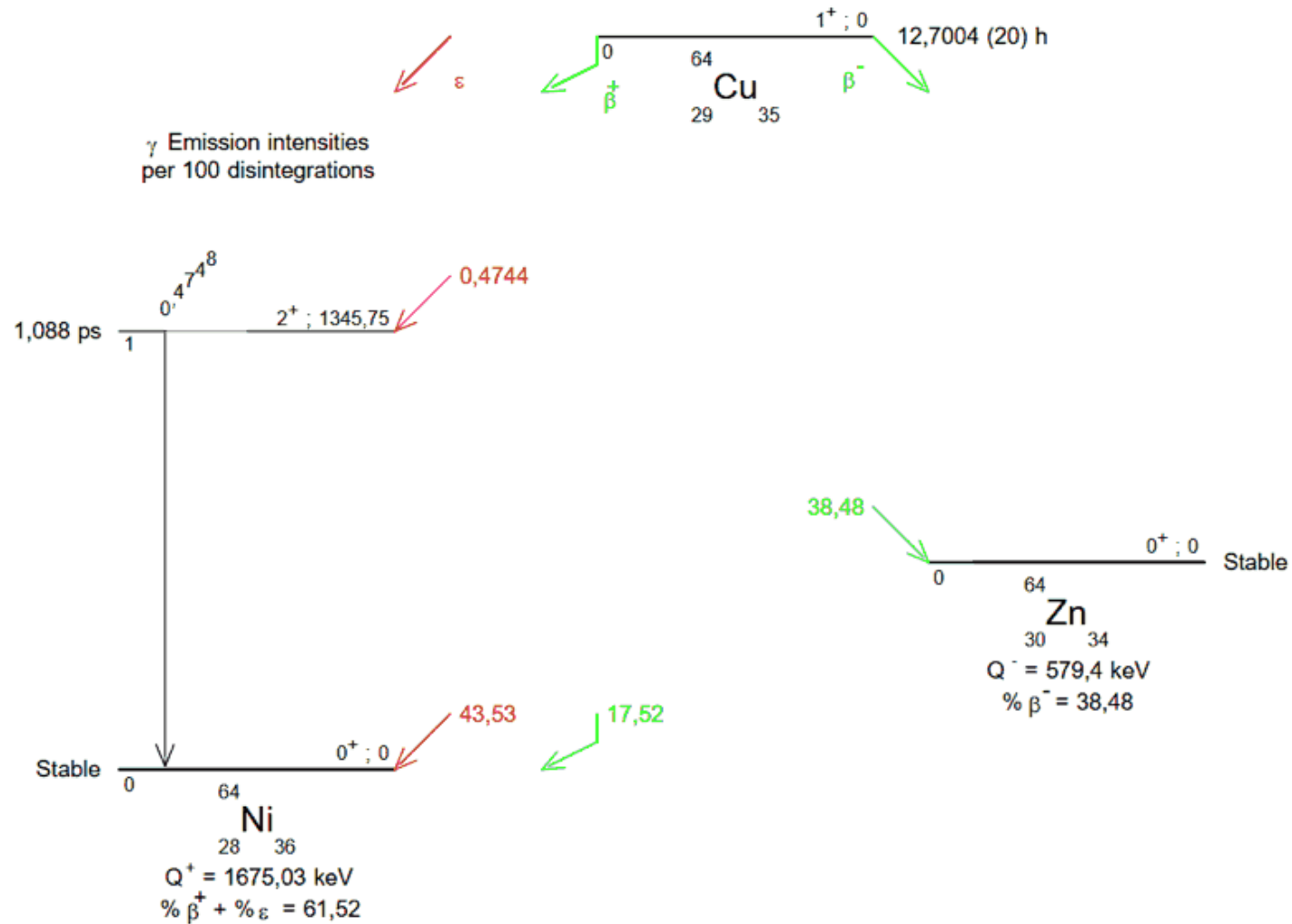

## EGS\_RadionuclideSource (cont.)

```
:start source:
  name          = my_mixture
  library        = egs_radionuclide_source

  base source    = name of the source used to generate decay locations

  activity       = [optional, default=1] total activity of mixture,
                  assumed constant. The activity only affects the
                  emission times assigned to particles.

  charge         = [optional] list including at least one of -1, 0, 1, 2
                  to include electrons, photons, positrons and alphas.
                  Filtering is applied to ALL emissions (including
                  relaxation particles).
                  Omit this option to include all charges - this is
                  recommended.

  experiment time = [optional, default=0] time length of the experiment,
                  set to 0 for no time limit. Source particles generated
                  after the experiment time are not transported.

:start spectrum:
  # See next slide...
:stop spectrum:
```



# EGS\_RadionuclideSpectrum

```
:start spectrum:
  type          = radionuclide
  nuclide        = name of the nuclide (e.g. Sr-90), used to look up the
                  ensdf file as $HEN_HOUSE/spectra/lnhb/ensdf/{nuclide}.ensdf
                  if ensdf file not provided below
  ensdf file     = [optional] path to a spectrum file in ensdf format,
                  including extension
  atomic relaxations = [optional, default=eadl] eadl, ensdf or off
                  By default, 'eadl' relaxations use the EGSnrc
                  algorithm for emission correlated with
                  disintegration events. Alternatively, 'ensdf'
                  relaxations statistically sample fluorescent
                  photons and Auger emission using comments
                  in the ensdf file. Turning this option off
                  disables all relaxations resulting from
                  radionuclide disintegration events.
  output beta spectra = [optional, default=no] yes or no
                  whether or not to output beta spectra to files.
                  Files will be named based on the nuclide and
                  maximum energy of the beta decay:
                  {nuclide}_{energy}.spec
# See next slide...
```

# EGS\_RadionuclideSpectrum (cont.)

alpha scoring = [optional, default=none] none or local  
Whether or not to deposit alpha particles locally.  
Since alpha particles are not transported in EGSnrc,  
there are only two options. Either discard the alpha  
particles and their energy completely, or deposit  
the energy immediately after creation in the  
local region.

extra transition approximation = [optional, default=on] on or off  
If the intensity away from a level in a radionuclide  
daughter is larger than the intensity feeding the  
level (e.g. decays to that level), then additional  
transitions away from that level will be sampled if  
this approximation is on.  
They will not be correlated with decays, but the  
spectrum will produce emission rates to match both  
the decay intensities and the internal transition  
intensities from the ensdf file.

:stop spectrum:  
:stop source:

# Source normalization

Sources provide a method for extracting the final fluence via the `getFluence()` method. However, the returned value depends on the source type:

- `EGS_IsotropicSource` returns the number of histories  $N$
- `EGS_PointSource` returns the number of histories  $N$
- `EGS_BeamSource` returns the number of histories  $N$
- `EGS_PhspSource` returns  $N_{\text{read}}/N_{\text{particles}} \times N_0$ , where  $N_0$  is the number of primary histories and  $N_{\text{particles}}$  number of particles in phsp file.
- `EGS_ParallelBeam` returns  $N/A$ , where  $A$  is the target shape area
- `EGS_CollimatedSource` returns  $N/d^2$ , where  $N$  is the number of particles hitting the target and  $d$  the distance defined by the user in the input file. If isotropic, you can set  $d = 1/\sqrt{2\pi}$  to recover a normalization to the number of histories. *Note that  $d$  is set to 1 by default!*
- `EGS_RadionuclideSource` returns the number of disintegration events

**It is important to be mindful of normalization when comparing results with other Monte Carlo simulation packages, or with experimental measurements!**

## ... more sources

- [EGS\\_SourceCollection](#) consists of arbitrary number of any other sources with user-defined weights.

```
:start source:  
    library = egs_source_collection  
    name = some_name  
    source names = list of names of previously defined sources  
    weights = list of weights for the sources  
:stop source:
```

- [EGS\\_AngularSpread](#) takes particle from other source and applies rotation to direction by angle sampled from a Gaussian distribution.

```
:start source:  
    library = egs_angular_spread_source  
    name = some_name  
    source name = the name of a previously defined source  
    sigma = positive: sigma in degrees  
            negative: fwhm in degrees  
:stop source:
```

**Finally, if you think this is not enough...**

⇒ implemented your own sources, deriving it from `EGS_BaseSource` or `EGS_BaseSimpleSource` class!

# **Browse the PIRS-898 html-documentation!**

<http://nrc-cnrc.github.io/EGSnrc/doc/pirs898/index.html>